



Domus 3

Dual therapy mode pressure redistribution system


Domus 3 combines alternating and static therapy modes to improve pressure relief care and reduce caregiver workload. It includes a low-pressure alert for patient safety and a weight indicator dial for easy adjustment. Additionally, it effectively wicks away moisture from the surface to aid in microclimate management, leading to a more comfortable experience and better prevention of pressure injury without losing cell firmness.

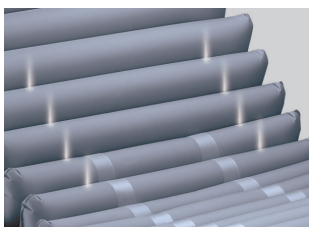

Microclimate Management
Micro low air loss effectively wicks away moisture from the surface, leading to a more comfortable experience.

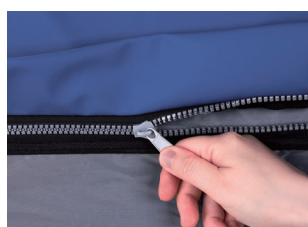

Decreases the Risks of Infection
The top cover is made with welded seams (high frequency induction hardening ®) to avoid fluid ingress.

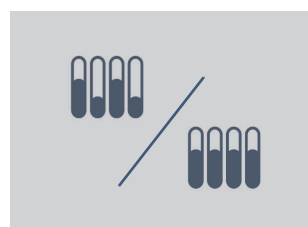

Dual Therapy Mode
Provide Static mode for enhanced patient immersion, Alternate mode for allowing pressure redistribution.



Prioritize Patient Safety
Visual signals are used to provide low-pressure alarms to alert caregivers in case of failure.

Domus 3

Revolutionize comfort, simplify care in pressure injury prevention



CPR Strap
Quickly deflate the mattress for emergencies.



Cable Management
Five loops on the side of mattress for managing cable to avoid tripping hazard.



Top Cover
Provide low shear, low friction, vapour permeable, and moisture protection with a basic level of flame retardancy.



Water Resistance Flame Resistance Biocompatibility Vapour Permeable



Quick Connector with Cap
Mattress can be transformed from a powered support surface to a sealed static air mattress for transportation need.

Specifications

Model	Domus 3				
Pump	Dimension	29 x 18.5 x 12.6 cm (11.4 x 7.3 x 5 in)			
	Weight	2.2 Kg (5.1 lb)			
	Case Material	Fire Retardant ABS			
	Supply voltage	120 V / 60 Hz			
	Therapy mode	Alternate / Static			
	Noise Level	< 35dB			
Mattress	Dimension	Cells 20	Length 200cm (78.7 in)	Width 90cm (35.4 in)	Height 20.3cm (8 in)
	Type	8" / Replacement			
	Weight	6.9 kg (15.2 lb)			
	Top cover material	Nylon / PU with Cotton Fabric			
	Cell material	Nylon / PVC			
	Maximum patient weight	200 kg (440 lb)			
	HCPCS Code	E0277			
	Ordering Information (For North America Only)				
	System	With Quilted Cover	PA01501		
Mattress	With Quilted Cover	PM09405			
Pump		PA01001			

Pump: water resistant standards (IPX0); Mattress: flame retardant standards (CAL117, EN597-1, EN597-2), RoHS, WEEE
Specifications are subject to change for improvement without notice

Wellell America Corp.
927 Mariner Street,
Brea, CA 92821, USA
T +1 714 671 3818
F +1 714 494 8173
info.usa@wellell.com
www.wellell.com
Print-2024/07 All rights reserved