



Procure Auto G3 Procure Auto Link G3

Maximize the efficiency of clinical workflow
with digitalized automation



Procure Auto G3 and Procure Auto Link G3 are advanced support surfaces designed to revolutionize patient care by leveraging the efficiency of clinical workflow.

The system combines intelligent pressure redistribution with automatic adjustments and cushion setting to deliver consistent, seamless support for patients at high to very high risk of pressure injuries.

Addressing inefficiencies in the support surface supply process



The duration of mattress preparation affects the waiting time, which increases the risk of pressure injuries.



The effectiveness of CPR is correlated with the time of emergency response.




Low utilization rate and maintenance efficiency in rental business.


Values and Solutions

Designed for professional healthcare facilities, addressing the needs of speed, reliability, and ease of use for caregivers.

Procure Auto G3 & Procure Auto Link G3 streamline the clinical workflow through digitalized automation, maximizing the operational efficiency to empower **Smarter, Safer, Faster** patient recovery care.

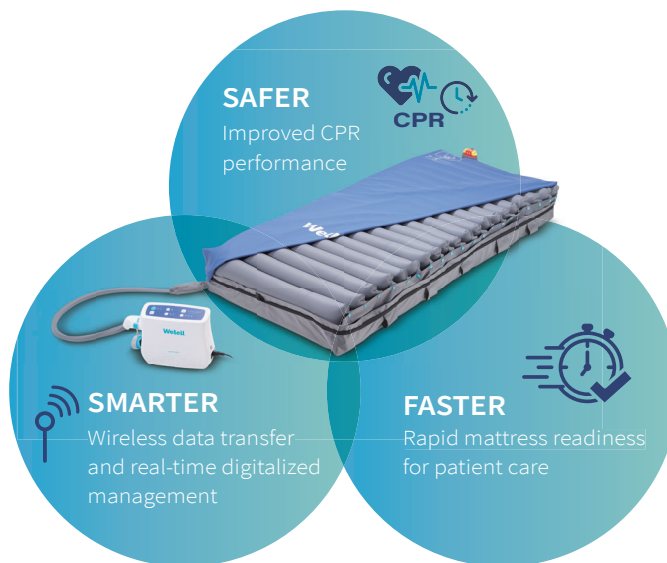

Link Connectivity

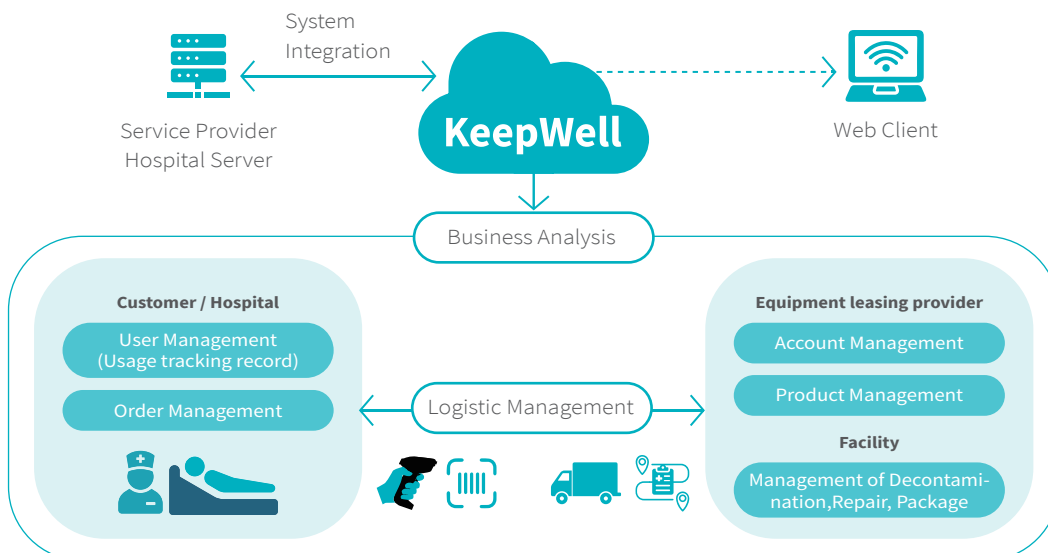


Procare Auto G3 & Procure Auto Link G3

Revolutionizing Patient Recovery Care with Smart Intuitive Design




Pressure Map

Subject: 181 cm/85 kg

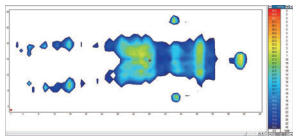

ALT pressure map

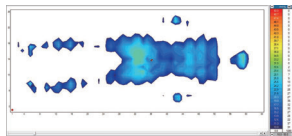

CLP pressure map

PUMP Features

- 1 Mattress & Cushion Integration**
The pump can be shared between the mattress and cushion, with automatic recognition for seamless system operation
- 2 Clear Icons for Intuitive Operation**
 - Max firm
 - Treatment modes : ALT / CLP(Static) / ALP (Sleep mode)
 - Adjustable comfort level / Audiovisual alarms
- 3 Auto Adjustment (AA)**
Automatically optimizes the pressure redistribution based on patient's weight
- 4 Wireless Data Transfer Function** 
Using 4G LTE communication for real-time tracking and data transmission
- 5 Twin-head Compressor**
Providing high flowrate to inflate the mattress within 11 minutes and operates at low noise.

MATTRESS Features

- 6 Press-Aid CPR Valve Enables Faster CPR Preparation**
An easy to access, press-to-release, single handed operation design performs quick deflation within 10 seconds
- 7 Easy Maintenance & Faster Assemble**
 - FlexVent (non-PVC air tunnels)
 - Color coded A/B snap-on tubes
 - Asset tracking tags
- 8 Durable Quick Connector**
 - A more durable fuel nozzle-style quick connector
 - Easy-to-clean rectangular air hose design

Europe Specifications

Pump



Procure Auto G3 & Procure Auto Link G3

Dimension (L x W x H)	12.7" x 10.4" x 5.2"			
Weight	7.9 lb			
Supply voltage	100V-120V, 60Hz			
Therapy mode	Alternating (ALT) / Continuous Low Pressure (CLP, Static) / Alternating Low Pressure (ALP) / Max Firm			
Cycle time	10 minutes			
Noise Level	< 22 dB **			
Communication (LINK)	Cellular : 4G-LTE, UMTS			
Dimension	Cells 20	Length 80"	Width 36"	Height 8"
Type	Replacement			
Weight	22lb			
Top cover material	4-way stretch PU, polyester cover with welded seams			
Cell material	TPU			
Maximum support weight	550 lb			

Mattress



Pump: water resistant standards (IP42); Mattress: flame retardant standards (EN597-1, EN597-2), WEEE

**Internal test conducted in a controlled environment

▲ Seat cushion is for add-on purchase