

Procare Elite

Alternating pressure redistribution system With low air loss


Procare Elite is suitable for patients at high to very high risk of pressure ulcers. This system offers ventilated low air loss for maximum comfort. Its cell-on-cell bladder design prevents bottoming out, while the unique heel relief function ensures comprehensive care.

Additionally, the system features five distinct treatment modes to meet the diverse needs of patients.



Alternating Pressure Relief

Pressure mapping reveals that in alternating mode, bladder A & B take turns to inflate and deflate can significantly relieve the body pressure.

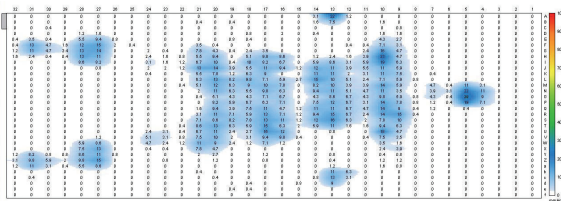


5 Modes to Fulfill Individual Needs

- Max Firm
- Alternating
- Continuous Low Pressure (CLP)
- Seat Inflation w/ Alternating
- Seat Inflation w/ CLP



4 Cycle time


10 Comfort Levels

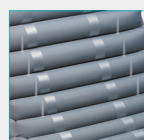

Ventilated
Low Air Loss
Wick moisture away to keep the patients' comfort




Cell-on-Cell Bladder
The bottom cells act as basic pressure support to avoid bottoming out in case of power failure

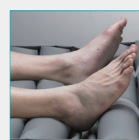


Heel Relief Function
Single air-cell deflation can float the heels with zero pressure



CPR Valve
Quickly deflate the mattress during emergency

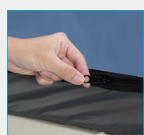

Zipper-all-around Cover
Reduce fluid ingress and simple one step to remove for cleaning

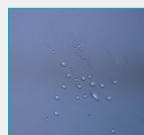

Quilted cover

- Low friction and shear
- Vapour permeable
- High-flame retardant
- Water resistant



Biocompatibility Vapour Permeable Water Resistance Flame Resistance



Cable Management
For a safe operational environment



Quick Connector
Once with cap on, instantly switch to static air mattress for transportation need



Visual and Audible Alarm

- Low Pressure
- Power Failure
- Tech Support



High Compressor Volume Airflow

Specifications

Model	Procure Elite	36"	42"	48"
Pump	Dimension	13.4" x 5.9" x 9.2" / 34.0 x 15.0 x 23.4 cm		
	Weight	6.83 lbs / 3.1 kg		
	Case Material	Flame retardant ABS		
	Supply Voltage	AC 110-120V 60 Hz, 0.17A		
	Operating Cycle	10/15/20/25 min. selectable		
Mattress		8' Mattress		
	Dimension	78.7" x 35.4" x 8"	78.7" x 42" x 8"	78.7" x 48" x 8"
		200 x 90 x 20.3 cm	200 x 107 x 20.3 cm	200 x 122 x 20.3 cm
	Cell Height	8" x 20 cells		
	Weight	16.5 lbs / 7.5 Kg	18.8 lbs / 8.5 Kg	21.0 lbs / 9.5 Kg
	Material	Cover: Non-Stretch Quilted cover; Cell: Nylon/TPU; Base: Nylon/PVC		
	Maximum Patient Weight	440 lb / 200 kg	1000 lbs / 455 kg	1000 lbs / 455 kg
	Flame Retardant Standards	CAL 116 & CAL 117		
	HCPCS Code	E0277		
Ordering Information (For North America Only)				
System	With Quilted Cover	PD04501	PD04504	PD04507
Mattress	With Quilted Cover	PM09407	PM09501	PM09601
	With Non-Quilted Cover	PM09408	PM09502	PM09602
	With Side Bolster, Quilted Cover	PM09409	PM09503	PM09603
Pump			PD04001	

Specifications are subject to change for improvement without notice

Wellell America Corp.

927 Mariner Street,
Brea, CA 92821, USA
T +1 714 671 3818
F +1 714 494 8173
sales.usa@wellell.com
www.wellell.com
Print-2024/05 All rights reserved